


FTIR Sampling Accessory for JASCO FT/IR-4000/6000 Series **DR PRO410-M** Diffuse reflectance accessory

The DR PRO410-M measures the diffuse reflectance of incident radiation focused on a sample. It is possible to measure 13 samples with a background by keeping purge. Applications include powder samples, rough solid surfaces, micro-samples, qualitative analysis of HPLC elutions and the study of substances adhered to powder surfaces.

14-position sample holder
Full purge capability
Coated mirror surface

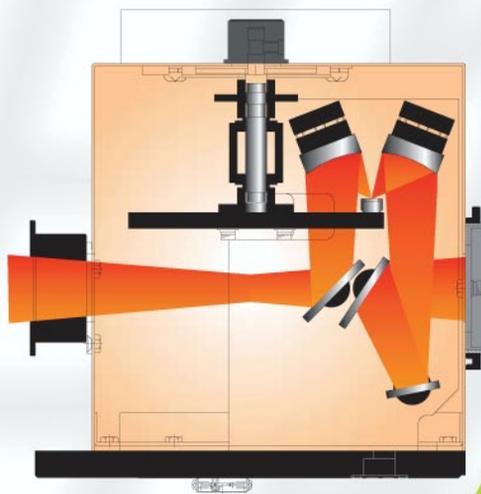

Applications

Materials which cannot be measured by the normal KBr method (e.g., clay, cement, minerals, pigments).
Analysis of minute samples
Qualitative analysis of HPLC fractions
Research on materials adhering to and absorbed by the surface of powders

FT/IR-4000 Series
Fourier Transform Infrared Spectrometer
IQ Accessory Recognition
Quick start external scan button



JASCO